


Dual Mono-Stereo Adapter

USER QUICK GUIDE

Introduction

The **Dual Mono-Stereo Adapter** enables two separate switches to be connected into a single stereo input. This is equivalent to a 3.5 mm TRS to Dual 3.5 mm TSF Stereo Breakout Cable.

Features

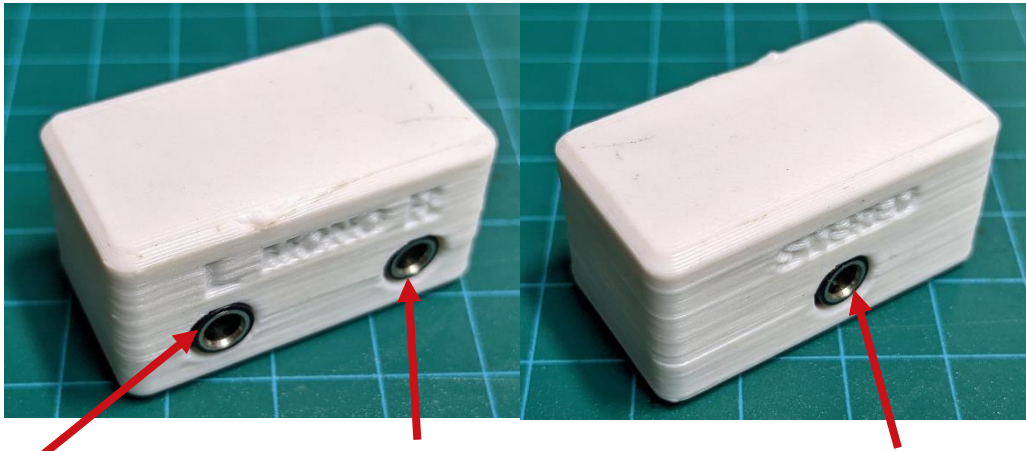

Left Mono Input

Right Mono Input

Single Stereo Output

Usage

Connect the switches to the desired left or right input jack. Use a 3.5 mm extension cable to connect the output to the input on the assistive device.

Compatibility

Switches: Assistive switches and devices with 3.5 mm output plug.

AT Devices: Assistive devices with 3.5 mm stereo input plug (e.g., Origin SWIFTY, Origin TAPIO, Origin HEADMOUSE, Darci USB)

Specifications

Mass [g]	12
Dimensions [mm]	38 X 21 X 19

Cleaning

The Dual Mono-Stereo Adapter can be wiped down with sanitizing cleaners.

Do not submerge in liquid.

Do not expose to high heat.